


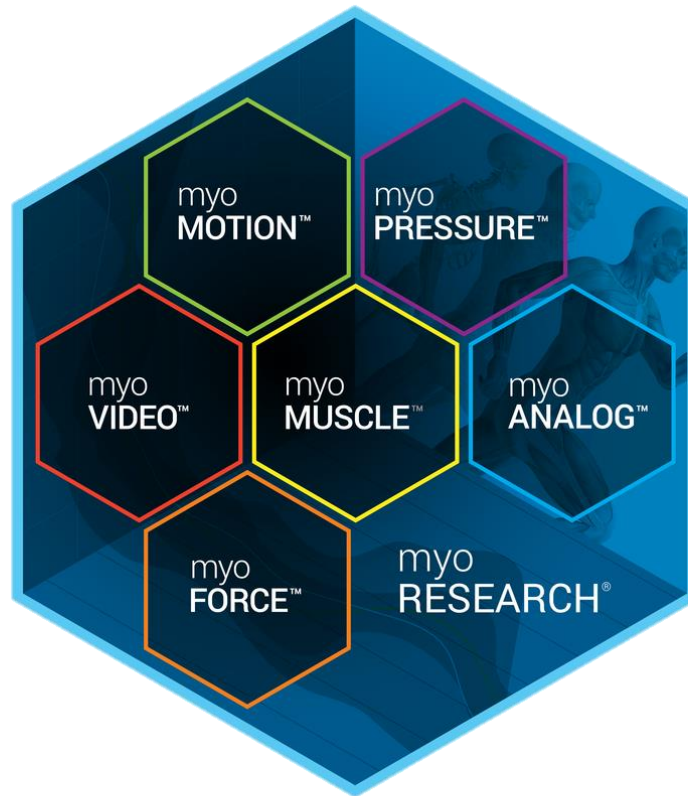

myoRESEARCH® 3.16

User Manual

Table of Contents

1	Introduction.....	5
1.1	Brief Description:	5
1.2	Intended Use	6
2	Computer Requirements	7
3	Software Security	8
3.1	Captured Data	8
3.2	User Entered Information	8
3.3	Stored Data	8
3.4	Data Storage Location.....	8
3.5	Methods to Access Data	9
3.6	Data Export Options	9
3.7	Safety Information Summary.....	9
3.8	System Level Options	9
4	Software Installation.....	10
4.1	Installation	10
4.2	Activation	10
4.3	Updates	10
5	Software Features Overview.....	11
5.1	General Software Navigation	11
5.2	Multi-device integration capabilities	12
5.3	Integrated Help System.....	12
5.4	Importing, Exporting, Streaming, and Backing up Data	12
6	Home	13
6.1	Hardware Setup.....	13
6.1.1	Insert and Configure Hardware	13
6.1.2	Save Hardware Profiles.....	14
6.2	Software Setup	14
6.2.1	About & License	14
6.2.2	Common	15
6.2.3	HTTP Streaming.....	15
6.2.4	myoMOTION Software Settings	15
6.3	Setup for a new measurement	16
6.4	Measurement Options Tool Bar	16
7	Measure	17
7.1	Record a Measurement.....	17
7.2	Actions Toolbar.....	17
7.2.1	Hardware-Related Measurement Tools	17

myoRESEARCH 3.16 User Manual

7.2.2	New Marker – Real Time Marker Menu	18
7.2.3	Screen Layout	18
8	Database	19
8.1	Database File Organization	19
8.2	Database List Functions	19
8.3	Selecting a Record to View	19
8.4	Database Tools	19
8.4.1	Database Type	19
8.4.2	Actions	20
9	Viewer	21
9.1	Main display features of the Viewer	21
9.2	Viewer Operations Toolbar	23
10	Report	24
11	Support	24
11.1	Submitting Technical Support Requests	24
11.2	Website Link to FAQ	24

- Manufacturer:

Noraxon U.S.A. Inc.
15770 North Greenway-Hayden Loop, Suite 100
Scottsdale, AZ 85260
Tel: (480) 443-3413
Fax: (480) 443-4327
Email: info@noraxon.com
Support Email: support@noraxon.com
Web Site: www.noraxon.com

- Authorized European Representative:

EC	REP
-----------	------------

Advena Limited, Tower Business Centre, 2nd Flr., Tower Street,
Swatar, BKR 4013 Malta

Website: <http://www.advenamedical.com>

CE Marking & Notified Body:

Clearance to market this product (#402, #405) in the European Community has been certified by
Notified Body #2797, British Standards Institution (BSI)

CE 2797

No part of this document may be copied, photographed, reproduced, translated, or reduced to any electronic medium or machine-readable form without prior written consent of Noraxon U.S.A. Inc.

Noraxon and myoRESEARCH are registered trademarks and the Noraxon logo, myoANALOG, myoFORCE, myoMETRICS, myoMOTION, myoMUSCLE, myoPRESSURE, myoVIDEO, myoSYNC, NiNOX, TRUsync and Ultium are common-law trademarks of Noraxon U.S.A., Inc. All other trademarks are the property of their respective owners. ©2018, all rights reserved. NiNOX and TRUsync are common-law trademarks of Noraxon USA in the USA and other countries. All other trademarks are the property of their respective owners. © 2020, all rights reserved.

1 Introduction

1.1 Brief Description:

The myoRESEARCH software suite (also known as MR3) fully integrates between the different software module options. Therefore, all software modules share the same basic menu operations. These operations are detailed within this manual.

For additional software tools, options, and suggestions for best practices, visit www.noraxon.com/courses.



1.2 Intended Use

Noraxon's myoRESEARCH software suite (MR3) is intended to compile, synchronize, and display biomechanical data that is measured by external, compatible, hardware sources. The benefit provided with use of MR3 is the provision of objective measures, from external hardware, that can assist the assessment of pathological human movement conditions and gauge any subsequent improvement offered by therapy, training or design changes. MR3 does not diagnose, but it does provide a platform for visualizing data obtained from external sources that are only to be used by qualified business professionals (see 'intended users').

Intended Users

Researchers or individuals trained in physical medicine, physical therapy, human performance or ergonomics and who have gained familiarity with this software and any applicable hardware.

Subject Populations – Medical

Individuals with neurological disorders, physical injuries, pre/post-surgical or post stroke conditions.

Subject Populations – Non-medical

Athletes, workers at their work site, subjects in new product trials.

Common Applications

Gait analysis; identification of inconsistencies and abnormalities; tracking over time the outcome of surgical, therapeutic or orthotic interventions; identification of ergonomic stress factors in the workplace or new product designs.

Suitable Environments

MR3 is suitable for use in hospitals, clinics, academic laboratories except near active HF SURGICAL EQUIPMENT or MRI IMAGING EQUIPMENT.

Expected (Essential) Performance

MR3 synchronously acquires and displays signals from one or multiple hardware source(s) in real time.

Consideration should be given to common place wireless networks that may be present in the environment. Under the conditions of excessive wireless communications performance of MR3 may exhibit sluggish (non-real-time) behavior.

Risk-Benefit

There is **no identified risk of physical harm or injury** with use of MR3. The benefit provided by use of the software is the provision of objective measures to assess the severity of pathological human movement conditions and gauge any subsequent improvement offered by therapy, training or design changes.

Precautions

Patients with epilepsy may experience an adverse response to the visual display (flashing lights / computer screen). Such patients should be treated with caution and careful monitoring during use.

2 Computer Requirements

The following computer system specifications are required for optimal use of Noraxon's myoRESEARCH software platform and modules and must be met for any multi-device MR3 installation. Computer systems sold and shipped by Noraxon will meet or exceed these specifications, and software compatibility is warranted for three months after the date of the Sales Quote. These specifications are subject to change without notice; please always check with your Noraxon salesperson for your particular use case.

O/S: Windows 10, 64-bit

Processor: Intel® Core™ i7 Processor (8th generation, PassMark score 10000 or better)

Memory: 16 GB RAM (minimum)

Primary Hard Drive: 256 GB SSD (M2 NVMe, larger capacity optional) ¹

Secondary Hard Drive: 1 TB SSD or 7200 rpm HDD (recommended for database)

Graphics Processor: A GPU with PassMark score 3000 or better is required (Nvidia Quadro P600 or GTX1050 GPU or better is recommended) ²

OpenGL: Version 4.5 or higher

OpenCL: Required

USB: A minimum of 2 USB 3.0 (Type A) ports ^{3, 4}

Camera: Noraxon NiNOX or Logitech C930e equivalent if the camera is to be used as a reference/tracking camera

When using NiNOX™ cameras, Noraxon USB cables must be used. Each pair of NiNOX Cameras require a separate USB 3.0 controller.

Portable Lab 2 Users: Computer must have 1, USB3 (Type-C) port and 1, USB3 (Type-A) port. An additional USB3 (Type-A) port will be required for multi-device systems (e.g. addition of force plates, pressure platform, etc.) ⁵

Recommended Computers for Multi-Device Systems:

Laptop (*up to 2 NiNOX™ cameras*)

- Dell Precision (i7 CPU, 16 GB RAM, 1TB M2 NVMe, NVidia Quadro P600)

Desktop (*up to 4 NiNOX™ cameras depending on multi-device configuration*)

- Dell XPS Tower (i7 CPU, 16 GB RAM, 1TB M2 NVMe, NVidia GTX1060)
- Dell Precision (i7 CPU, 16 GB RAM, 1TB M2 NVMe, NVidia Quadro P600)

¹ Noraxon recommends use of an external hard drive for data backup.

² Noraxon highly recommends that you update your graphic board drivers using the manufacturer's website.

³ One port is needed for each Noraxon device to be connected. A Noraxon provided USB hub can be used to add ports. More than one USB controller may be required when combining myoVIDEO and myoPRESSURE.

⁴ Note the cable length limit is 5 meters for USB 2.0 and USB 3.0. Longer lengths can be accommodated using Noraxon provided active extension cables or USB hubs.

⁵ Refer to the Portable Lab Quick Start Guide for additional settings for Portable Lab or Portable Lab 2 system.

3 Software Security

3.1 Captured Data

Noraxon's MR3 software is used to capture and display video, noninvasive electrophysiological, kinetic, and kinematic signals. All Noraxon hardware devices communicate directly to the collection computer through standard USB 2.0 and USB 3.0 connections. No form of encryption is used for data transfer as this data is not considered sensitive under HIPAA rules and HITECH regulations (see list of PHI data elements below).

- Names
 - Telephone numbers
 - Fax numbers
 - Email addresses
 - Social Security numbers
 - Medical record numbers
 - Health plan beneficiary numbers
 - License plate numbers
 - URLs
 - Full-face photographic images
 - Any other unique identifying number, characteristic, code, or combination that allows identification of an individual
- Please note that video capture is the only measurement modality that can result in sensitive patient information. It is the responsibility of the system users to not record the patients face when conducting kinematic video measurements.

3.2 User Entered Information

All users that have default level access to the operating system running Noraxon's MR3 software will also have permissions to write to the Noraxon MR3 software folder.

Each Subject entered in the software has information fields available to enter data that would be considered PHI data elements, such as name and address. Patient information considered to be sensitive should not be entered into Subject demographic. It is the software user's responsibility to avoid these information placeholders when HIPAA rules are applicable.

To better suit the needs of various Noraxon customers, the software data fields are completely customizable and can be restricted to satisfy the needs of the institution.

3.3 Stored Data

Subject demographic data, as determined by institution needs. Acquired data, both as measured variables and video streams.

Definitions of acquired signals taken in the measurements include electromyography, kinetic, video and inertial sensor based kinematic data. All data recorded, in and of itself, is not a means to diagnose care but must be interpreted by a professional with access to patient medical records.

3.4 Data Storage Location

Data collected using the Noraxon hardware is stored locally on the machine running Noraxon's MR3 software.

myoRESEARCH 3.16 User Manual

An external database location can be configured within the MR3 software (MR 3.10 or higher), but the user should understand the limitations of data transfer (i.e. transfer speed) depending on the designated location.

The Noraxon software database consists of three tiers including the Project Name, Subject Name and Record Name. The data is not encrypted.

3.5 Methods to Access Data

The primary MR3 software (used to collect data directly) can be used to: display the data graphically in real-time and post-collection, generate summary reports, and create video screen captures. Reports can be generated in PDF format and will contain only user-entered information.

3.6 Data Export Options

Raw data may be exported in various formats including Excel, MATLAB, ASCII (CSV), C3D or using the proprietary Noraxon database structure.

Using the Noraxon native format gives access to a graphical user interface for easy selection of records and reports.

3.7 Safety Information Summary

All data, including user entered PHI, is stored in a proprietary database accessible via Noraxon's MR3 software. The user should establish the necessary controls to restrict access to the MR3 software on the installed PCs.

Please note that nearly all above features can be uniquely removed from your installation image. The ability to record, save, export can all be tailored to your security needs. For example, you may like only to export using the Noraxon format and remove video reporting altogether.


3.8 System Level Options

The Noraxon MR3 software does not offer a direct link to an electronic patient record. Reports generated by the MR3 software can be in an electronic file format (e.g. PDF) for inclusion within EPR as a separate document.

4 Software Installation

4.1 Installation



1. Insert the MR3 USB flash drive into the PC
2. A menu will automatically pop up
3. Click on the Noraxon installation file and follow the Wizard's instructions

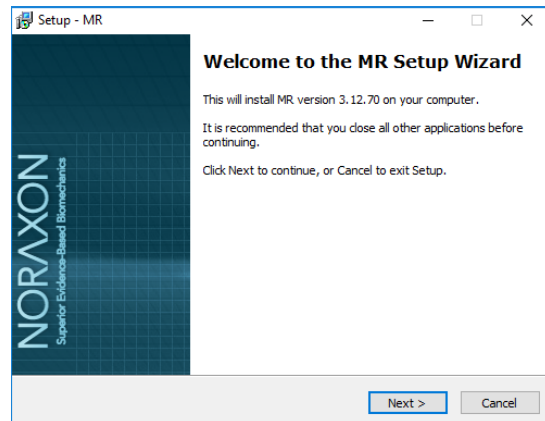

4.2 Activation

The installed companion software must be *activated* before unrestricted use is possible.

1. Open MR3
2. A dialog box will indicate how many more times MR3 can be opened
3. Click on "Activate"
4. Enter the License ID provided on your USB flash drive and press "OK"
5. If you have an internet connection, click Activate by Internet for immediate activation
6. Alternatively, email the provided activation ID to activation@noraxon.com. Noraxon Support will email or respond by phone with the Activation Code. Enter the provided Activation Code to remove any restrictions on use.

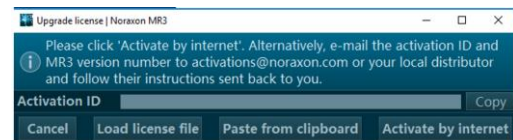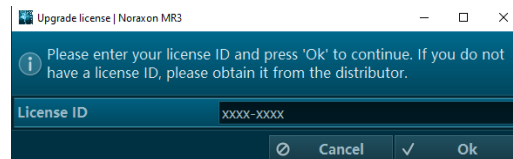


4.3 Updates

1. Perform a backup of the data folders to a separate drive as a precaution.
2. Click on the patch/update link provided via email or at noraxon.com/downloads.
3. Download the patch/update file.
4. To install the patch/update, click "Run" on the dialog box. No password is required.

5 Software Features Overview

5.1 General Software Navigation

Three major interface components are available on all screens.

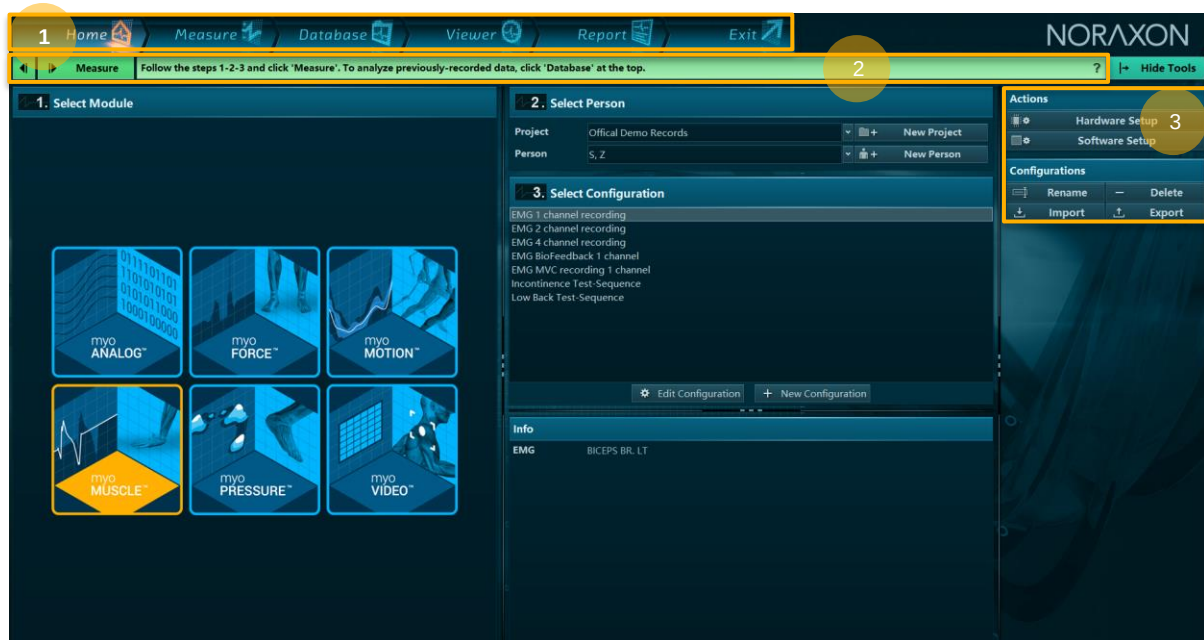

1. Main software navigation bar

Move between the stages of data collection by clicking through the options.

2. Green instruction line with main action button

Follow the instructions to operate to the next step in the workflow

3. Tool bar

The available options or settings for each screen are displayed here.


Select **Hide Tools** or **Show Tools** to display or retract this tool bar.

5.2 Multi-device integration capabilities

Your purchase of MR3 includes the software modules that correspond with your hardware purchase. The software modules you have access to will display on the Home screen. You will need to configure the hardware you intend to use in MR3 by following the Hardware Setup process described in the corresponding Hardware User Manuals. Once each device is configured in Hardware Setup, you can use it within a given measurement configuration. All devices will automatically be time synchronized with use of the Noraxon myoSYNC device.

5.3 Integrated Help System

Context specific help information is available by clicking the question mark.

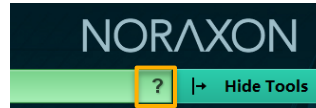

5.4 Importing, Exporting, Streaming, and Backing up Data

The basics of importing and exporting records with MR3 can be found in section 8.4 of this User Manual. For more detailed information and a step-by-step guide to perform each of these functions, please see the MR3 Import-Export feature guide.

6 Home

Before the MR3 software can record and store measurements, you must configure the installed software through the Hardware and Software setup and prepare the recording configurations for the measurement. These functions are available from the Home screen.

6.1 Hardware Setup

6.1.1 Insert and Configure Hardware


This step is only needed when new devices are being added to the software or a change to the hardware profile is needed.

Step 1

Connect all measurement devices to the PC.

Step 2

Select **Hardware Setup**

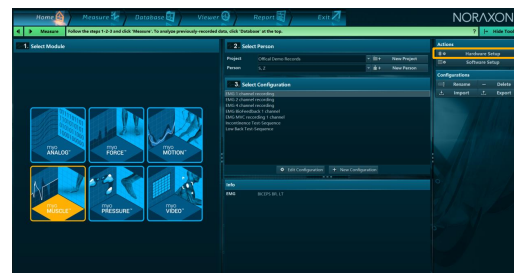


Step 3

Drag and drop (or select **Insert**) to move the hardware of interest to the Selected devices window from the New devices window.

Step 4

Select **Configure** to access device specific settings and individual sensor detection.

Each device has a specific Configuration menu. Consult the corresponding Hardware User Manual for the specific hardware configuration procedures for your device.

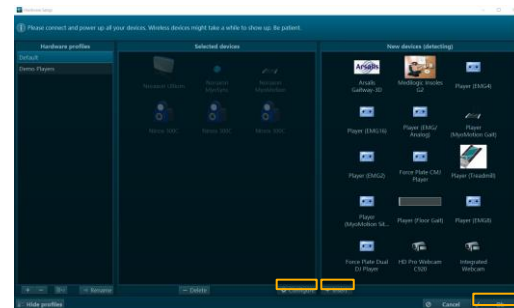

Step 5


Select **Ok** to exit Hardware Setup.

6.1.2 Save Hardware Profiles

Multiple hardware profiles can be saved. This allows specific hardware setups and device configurations to be saved for later use. This is very useful if more than one EMG and biomechanical sensors are combined in certain multi-device setups.

Step 1

Select **Profiles** within the Hardware Setup menu to view saved hardware configurations.

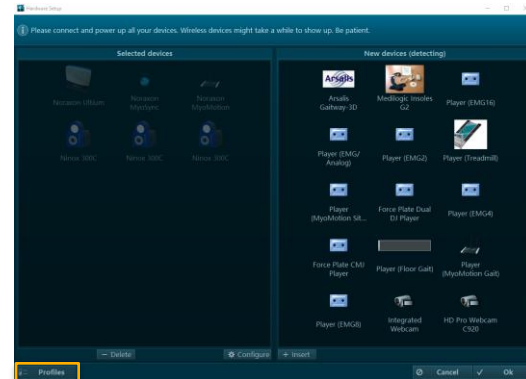


Step 2

To create a new profile, copy, delete, or rename an existing profile select the corresponding tool at the bottom of the Hardware profile window.

Functions from left to right:

- Create a new profile
- Delete marked profile(s)
- Copy selected profile
- Rename selected profile

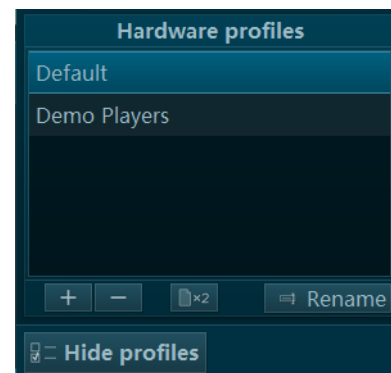


6.2 Software Setup

In the Software setup you can access general software settings as well as some module specific functions. The basic settings needed to operate the software are described here.

6.2.1 About & License

Information about the software version, license, and manufacturer can be found in the **About & License** section of the Program Settings window.

If you wish to upgrade your existing license with a new module (e.g. if a myoMUSCLE-only license needs to be upgraded to also include the myoVIDEO module), click **Upgrade** and follow the on-screen instructions. Select **Deactivate** to remove the access rights from the PC.

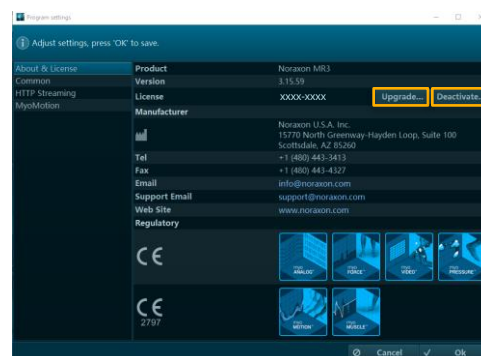

6.2.2 Common

Basic software settings can be adjusted in the **Common** section.


Select the operating **Language** using the dropdown menu (default language is the language of installed Windows OS). Languages supported by MR3 are: English, German, Chinese, and Japanese.

Select the **Screen font DPI** using the dropdown menu and data field.

Establish an **External database folder** by selecting the checkbox.


Select video compression settings

Set a password to protect access to the MR3 software

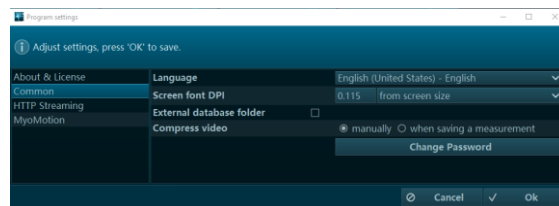


6.2.3 HTTP Streaming

Enable real-time data streaming for 3rd party programs and set the access port in the **HTTP Streaming** section.

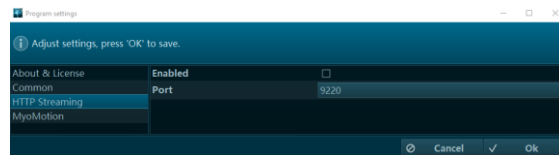

6.2.4 myoMOTION Software Settings

This section is only operational for the myoMOTION 3D inertial sensor system. Ignore it in case of using other systems without the myoMOTION sensors. See the myoMOTION User Manual for specifics.

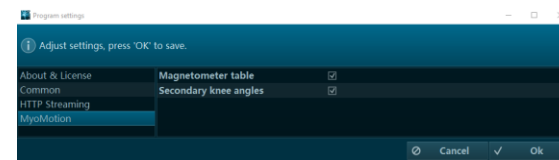

6.3 Setup for a new measurement


Step 1

Within the Home screen, click on the software module icon that corresponds with the device(s) in your measurement.

Create a **New Project** or select an existing project.

Create a **New Subject** or select an existing subject.

Select the configuration you would like to measure within. If you need to make changes or view the devices within a configuration select **Edit Configuration**. Select **New Configuration** to create a new measurement configuration.

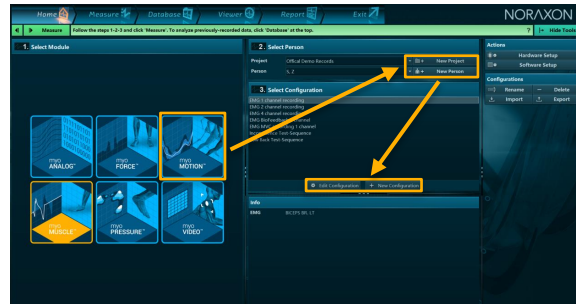


Step 2

Insert the devices to be used for the measurement into the configuration by dragging a device in from the list of **Available Devices**.

Use the sensor list to select or deselect the number of sensors and/or measurement channels to be used within the configuration.

Each system has a variety of unique measurement settings and option. Refer to device specific Hardware User Manual for the details for each system.

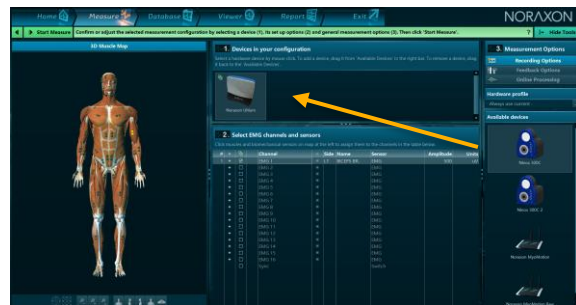

Step 3

Save the configuration.

To save your configuration you have two options:

- Select **Home**. This will return you to the Home screen and prompt you to save or discard the configuration changes.
- Select **Start Measure**. This will move you to the Measurement screen and save your configuration changes.




6.4 Measurement Options Tool Bar

Within the Edit Configuration window, you can edit advanced measurement options.

myoRESEARCH 3.16 User Manual

To learn more about operating one of these advanced settings, select the integrated help icon or visit www.noraxon.com/courses.

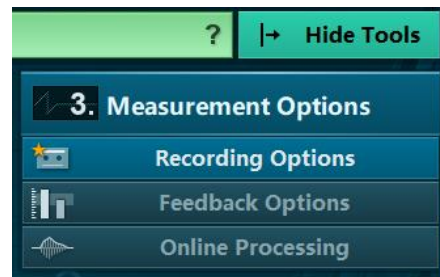

7 Measure

7.1 Record a Measurement

To record a measurement or view your data in live mode, make sure you have first set up your measurement configuration as outlined in Section 6 of this manual.

Step 1

Following subject and measurement configuration setup procedures, select **Measure** to move to the measurement screen.



Prior to starting a recording, determine if your measurement devices require specific calibration steps. See the device's corresponding Hardware User Manual for an overview. If required, perform the appropriate calibration procedures.

Step 2


Check the signals from the sensors and, if acceptable, follow the steps in the green tool bar.



Select the back arrow (◀) to return to the **Edit Configuration** window and adjust the configuration settings and measurement options.

Step 3

Upon stopping a recording, save or discard the measurement using **Discard & measure again**, **Save & measure again**, or **Save & view**

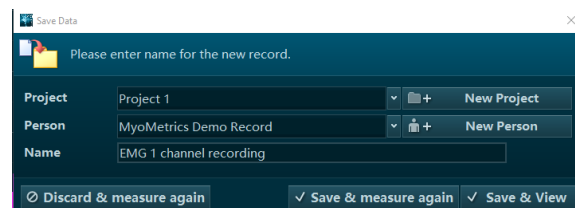

7.2 Actions Toolbar

7.2.1 Hardware-Related Measurement Tools


Depending on the hardware you are using you will see a variety of unique measurement tools. Make sure to reference your device Hardware User Manual for more information on operating each tool.

7.2.2 New Marker – Real Time Marker Menu

With this menu you can create marker names and manually mark or comment on events in your recording.

Step 1


Enter a marker name and press **Add** to add it to the list of available markers

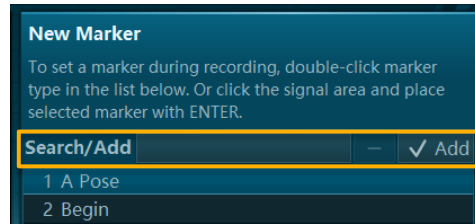

Step 2


When recording, the **Add** button changes to **Set Marker** button.

To place a custom marker type in your recording, first click on the desired marker type in the list then click on the **Set Marker** button in marker menu. An acoustic beep will confirm when a marker is placed in the recording. Alternatively, the keyboard **Enter** button can be pressed to set a marker.

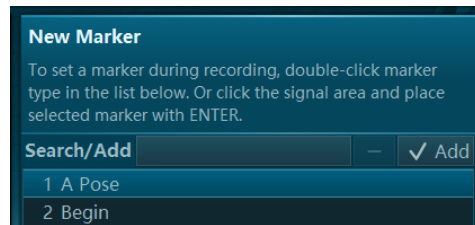

Additional Option

Import, Export, and Delete buttons allow you to edit or import/export the list of markers (from or to other MR3 versions).

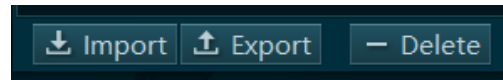


7.2.3 Screen Layout

This list of tools will allow you to customize the format in which your data is displayed during measurement and review.

The grouping options provide default reorganization of the signals displayed in each window.

The interaction commands (e.g. copy, paste, save, etc) allow for custom layouts to be used and shared as needed.

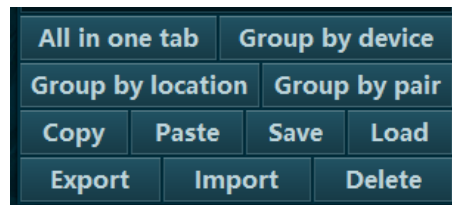


8 Database

The Database screen manages all functions related to data file management such as selection, editing, exporting, and importing.

8.1 Database File Organization

The Database file system is organized in a hierarchical system as follows.


1. Projects
 - a. Subjects
 - i. Records

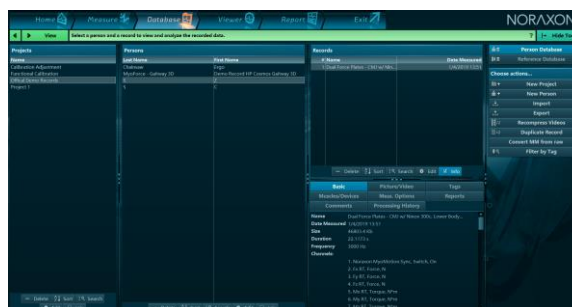

8.2 Database List Functions


At the bottom of each column of the database there are a series of functions that allow for typical operating procedures within the Database, such as:

- **Delete** (Remove project, subject, or record from the Database)
- **Sort** (Sort Database in ascending or descending order)
- **Search** (Quickly locate a project, subject, or record in a large database)
- **Edit** (Change or add information stored at each level of the database)
- **Info** (Use to view available information on a project, subject, or record)

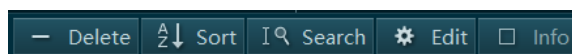

8.3 Selecting a Record to View

To select a record to view either double-click on the record of interest or select **View**.

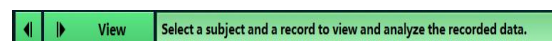


Two records can be viewed side-by-side by holding CTRL to select the two records of interest and clicking **View**.

8.4 Database Tools

8.4.1 Database Type


Subject Database

When Subject Database is selected you will see all projects, subjects, and records saved locally.

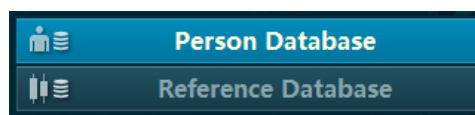

Reference Database

This database contains all reference data records that were compiled via the **Add to Reference** function in the Report screen.

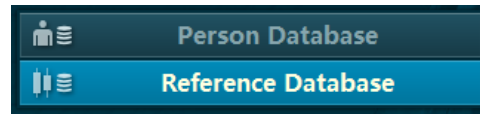

The reference data sample files included in your purchase of MR3 are not meant to be used for any medical diagnosis or similar purposes. It only serves as an educational example on how to use or operate normative data.

8.4.2 Actions

Several additional actions are available in the Tool window of the Database screen. See a description of each below.

- **New Project**
 - Create a new project in your database. Similar process to creating a new project from the Home screen.
- **New Subject**
 - Create a new subject in your database. Similar process to creating a new subject from the Home screen.
- **Import**
 - Import records from external locations such as previous MR software version or backup data from external hard drives or network folders.
- **Export**
 - Export data in native format to an **external location** on your PC or to another **Subject or Project** in your same MR3 database.
 - Export data in a different format; e.g. CSV file
- **Recompress Videos**
 - Reduces file size of collected raw video.
- **Duplicate Record**
 - Produces a replica of the record.
- **Filer by Tag**
 - Allows for quick searching for pre-tagged records.

When transferring data from one instance of MR3 to another, it is necessary to use the import/export "to/from external location" to achieve proper format for the data transfer or backup. For more detailed information and a step-by-step guide to perform each of the import/export functions, please see the MR3 Import-Export feature guide.

9 Viewer

Data loaded into the record Viewer screen will display in the same channel layout scheme that was used for recording. Within the Viewer screen you can inspect the quality and details of a record, apply signal processing, place markers for analysis period definition, and select additional Play Mode settings.

9.1 Main display features of the Viewer

Splitter Lines

Use the splitter lines to increase or decrease the size of the windows displayed in the Viewer to accommodate your preferences.

Device Tabs

For multi-device recordings, toggle between the tabs on the left side window to view available channels by device or custom configuration.




Right-click on a tab to access additional settings and operations.

Video/Animation windows

Video recordings and devices with animations (e.g. 3D motion, pressure, force plate CoP) will display playback visuals in the Video/Animation window.

Timeline & Playback Bar

The record timeline and playback options are located at the bottom of the Viewer screen.

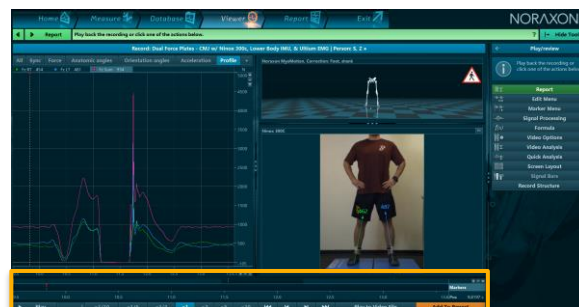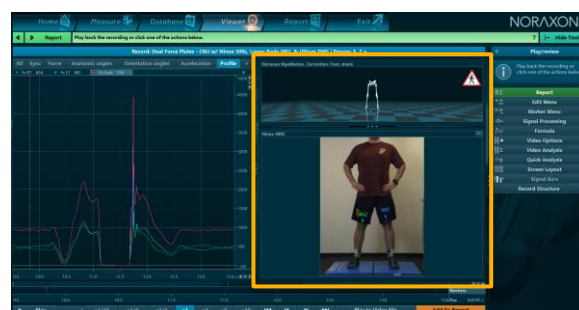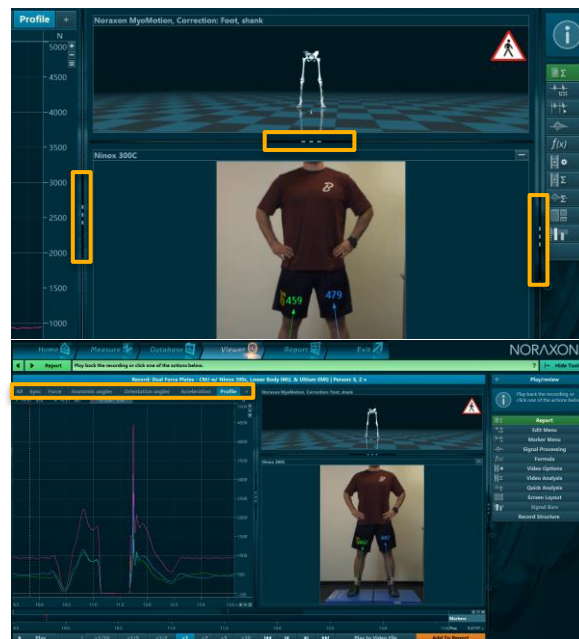

myoRESEARCH 3.16 User Manual

Signal Plot Display & Settings


Double-click or click-hold-and-drag the header of a signal to change its location.

Right-click on the header of a signal to access additional settings and operations.

Double-clicking the header will expand the data into its own window. Exiting the window will return the signal(s) to the original location.

Click, hold, and drag the header to place signal plots in different locations, orders, or tabs.

1. Bottom right arrow – places the dragged signal below the current signal plot in view
2. Top left arrow – places the dragged signal above the current signal plot in view
3. Centered arrows – overlays the dragged signal onto the current signal plot in view
4. It is also possible to drag and drop entire tabs of data following the same indicators above

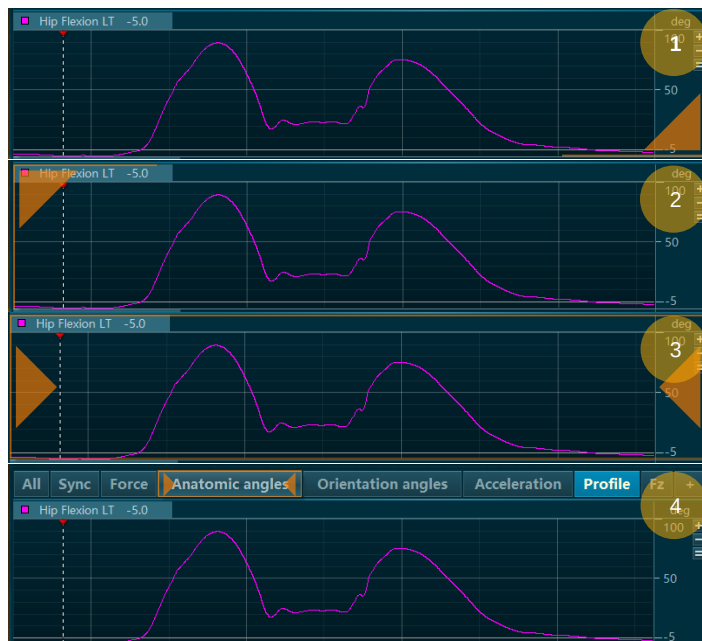

9.2 Viewer Operations Toolbar

The right-side toolbar is meant to make editing data in the viewer more user-friendly and accessible. It changes based on selections in the Viewer screen to display the corresponding features for the operation selected in the Viewer.

This is the initial view of the toolbar in the Viewer. Click **Show Tools** in the top right if the toolbar does not appear.

Report – executes reporting process; same action as Report button on green toolbar at the top of the screen

Edit Menu – allows user to modify and interact with recorded data

Marker Menu – allows user to edit and add markers to the recording

Signal Processing – opens signal processing menu

Formula – opens Python-based formula editor

Video Options – allows the user to modify view settings and interact with video recordings

Video Analysis – opens video analysis drawing tools; requires myoVIDEO module

Quick Analysis – executes a quick statistical report of recorded data for a selected portion of the recording. Includes a visual of the selected channel(s) and editable calculated parameters (e.g. mean, peak, minimum, area, etc)

Screen Layout – opens the screen layout options menu

Signal Bars – enables bar graph visualization of recorded data



10 Report

Several pre-installed reports are available for use within the Report screen. To access the Report screen, first open a record from the Database.

Reports are sorted in eight major tab sections. The accessibility of these tabs is dependent upon the equipment and modules installed.

- **myoMUSCLE** (Clinical, Master) – for EMG and other biomechanical sensor records
- **myoMOTION** – for IMU-based 3D kinematic records
- **myoPRESSURE** – for pressure distribution and gait analysis records
- **myoVIDEO** – for video analysis records
- **myoFORCE** – for 3D force and isokinetic analysis records
- **myoANALOG** – for analog records
- **Quick Analysis** – default report for the Viewer-based Quick Analysis

Report formatting can be edited to fit the needs of the user. For additional details on customizing reports, visit www.noraxon.com/courses.

11 Support

11.1 Submitting Technical Support Requests

A Support Request can be submitted using the online form available at this link:

<https://www.noraxon.com/support-learn/support-request/>

Provide all information requested by the form including a **detailed** description of the problem being experienced and your telephone number or e-mail address.

11.2 Website Link to FAQ

Answers to frequent questions can be found at Noraxon's Frequently Asked Questions (FAQ) website page at this link:

<https://www.noraxon.com/support-learn/technical-support/faqs/>

Other educational material is available at this link:

<https://www.noraxon.com/support-learn/technical-support/>